



**Fisherbrand**<sup>®</sup>  
QUALITY. RELIABILITY. VALUE.

# Liquid Handling Consumables


*Elite Pipettes and SureOne™ Tips*



This pipetting system provides a complete line of pipettes and tips with optimal fit. Take a look at the benefits that make the Fisherbrand™ Elite™ Pipette and Fisherbrand™ SureOne™ your high quality solution.

**New!**  
**SureOne Filter Tip Reload System**

Introducing the Fisherbrand Filter Tip Reload System that offers the same benefits as the SureOne pipet tips and is the only stackable filter reload system, providing an environmental solution without compromising sterility. With offset wells allowing the filter tips to stack, this advanced system reduces waste and saves bench top and storage space.

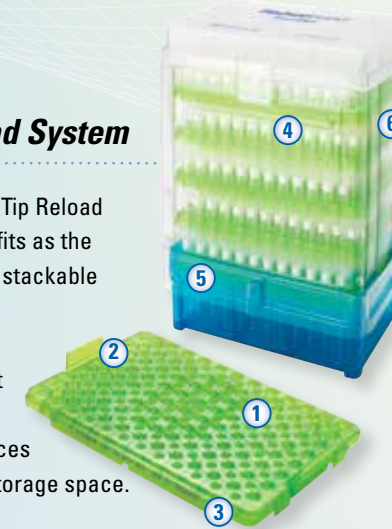

- 1 Offset wells allow filter tips to nest and stack.
- 2 Rigid inserts are multichannel format offering a robust non-pliable design providing secure tip seal.
- 3 Lot numbered inserts for traceability also include easy to read filter tip volume and Fisher Cat. No.
- 4 Sterile packaging with outer sleeve protects tips from exterior contamination while refilling inserts.
- 5 Environmentally-friendly filter tip packaging option reduces plastic waste and saves bench top space.
- 6 Blister sleeve guide rails align tip inserts for smooth transfer into the rack.

**Fisherbrand™ Elite™ Pipettes**

Achieve the best results in your lab with exceptional ergonomics while excelling at accuracy and precision. This comprehensive product line offers a wide range of pipettes for multiple applications providing an optimal fit with the SureOne pipet tips.



- 1 Superior Finger Rest
- 2 Index Finger Support
- 3 Hand Support
- 4 Optimized fit for SureOne Tips
- 5 Soft-touch Tip Ejector
- 6 Extremely Large Display
- 7 Large Surface Tip Ejector
- 8 Double Action Pipetting Button

**SureOne Rack**



- 1 Flush fit inserts sit flush with rack, allowing easy multichannel pipetting.
- 2 Recycle logo molded into all rack components, allows recycling companies to properly recycle the plastic.
- 3 Notched clasp designed to prop lid open for autoclave venting.
- 4 Sure stack lid rim prevents racks from sliding around when stacked and carted around a lab.
- 5 Safe lock latch allows for one-handed opening and closing ability and won't pop open when dropped.
- 6 Sure snap design securely holds insert in place, preventing it from being lifted out while pipetting.
- 7 Flip-lid can be either a hinged or lift-off lid; hinged format does not prevent use of a multichannel pipet.

**SureOne Reload Packaging**



- 1 Product information displayed on box.
- 2 Dust flap prevents debris from contaminating the product.
- 3 Finger knockouts allow easy access to slide out product.
- 4 Canyon cut package allows for easy tip insert removal.

**SureOne Standard Sterility Claims** Non-sterile products are certified to be free of RNase, DNase, DNA and Pyrogen. E-beam sterile products are certified to be free of RNase, DNase, Pyrogen, Bioburden, PCR Inhibitors, Endotoxin and DNA.

To order a SureOne Standards Certificate, contact Fisher Scientific Customer Service with the Cat. No. and Lot No.

# Fisherbrand® Elite Pipettes



The Fisherbrand™ Elite™ pipette has a new, improved ergonomic design providing index finger and hand support for maximum comfort. The low plunger force combined with a larger soft-touch ejection tip helps you work longer without fatigue.

Ensuring the success of the application where sterile products are required, the entire pipette can be steam autoclaved. The Elite pipette's double action pipetting button's rotating upper disc prevents accidental volume alteration during pipetting. An extremely large display provides better visibility of the selected volume. This comprehensive product line offers a wide volume range for multiple pipetting applications and excels in accuracy and precision.

- Soft-touch tip ejector
- Optimal fit to SureOne pipet tips
- Fully autoclavable
- Comfortable, lightweight handle with finger rest
- Extremely low plunger forces

## Variable Volume Single Channel Pipettes

Fisher Scientific Cat. No.	Description	Range	Increment	Accuracy	Precision	SureOne Tip for calibration
11815762	Fisherbrand Elite, 0.2-2 µl	0.2-2 µl	0.002 µl	2.5 to 12.0%	2.0 to 10%	11953416
11825762	Fisherbrand Elite, 0.5-5 µl	0.5-5 µl	0.01 µl	1.5 to 6.0%	1.0 to 5.0%	11953416
11835762	Fisherbrand Elite, 1-10 µl	1-10 µl	0.02 µl	1.0 to 2.5%	0.5 to 2.0%	11953416
11845762	Fisherbrand Elite, 2-20 µl	2-20 µl	0.02 µl	1.0 to 3.0%	0.4 to 1.5%	10678325
11855762	Fisherbrand Elite, 5-50 µl	5-50 µl	0.1 µl	0.6 to 3.0%	0.3 to 2.5%	10678325
11865762	Fisherbrand Elite, 10-100 µl	10-100 µl	0.2 µl	0.8 to 3.0%	0.2 to 1.0%	10678325
11875762	Fisherbrand Elite, 20-200 µl	20-200 µl	0.2 µl	0.6 to 1.8%	0.2 to 0.7%	10678325
11885762	Fisherbrand Elite, 100-1000 µl	100-1000 µl	1 µl	0.6 to 1.0%	0.2 to 0.6%	11568442
11895762	Fisherbrand Elite, 0.5-5 ml	0.5-5 ml	0.01 ml	0.5 to 2.0%	0.2 to 0.8%	11937754
11805772	Fisherbrand Elite, 1-10 ml	1-10 ml	0.02 ml	0.5 to 2.0%	0.2 to 0.8%	11957754
15268638	Fisherbrand Elite Kit: Includes four pipette sizes (0.2-2 µL, 2-20 µL, 20-200 µL, 100-1000 µL) and stand					



## Variable Volume Multichannel Pipettes

Fisher Scientific Cat. No.	Description	Range	Increment	Accuracy	Precision	SureOne Tip for calibration
11815772	Fisherbrand Elite, 8-Channel, 1.0 -10 µl	1.0-10 µl	0.02 µl	2.4 to 12%	1.6 to 8.0%	11953416
11835772	Fisherbrand Elite, 8-Channel, 5-50 µl	5-50 µl	0.1 µl	1.5 to 5.0%	0.7 to 2.0%	10678325
11825772	Fisherbrand Elite, 8-Channel, 10 -100 µl	10-100 µl	0.2 µl	1.3 to 5.0%	0.5 to 2.0%	10678325
11845772	Fisherbrand Elite, 8-Channel, 30-300 µl	30-300 µl	1.0 µl	1.0 to 4.6%	0.3 to 1.5%	11903456
11855772	Fisherbrand Elite, 12-Channel, 1.0 -10 µl	1.0-10 µl	0.02 µl	2.4 to 12%	1.6 to 8.0%	11953416
11875772	Fisherbrand Elite, 12-Channel, 5-50 µl	5-50 µl	0.1 µl	1.5 to 5.0%	0.7 to 2.0%	10678325
11865772	Fisherbrand Elite, 12-Channel, 10 -100 µl	10-100 µl	0.2 µl	1.3 to 5.0%	0.5 to 2.0%	10678325
11885772	Fisherbrand Elite, 12-Channel, 30-300 µl	30-300 µl	1.0 µl	1.0 to 4.6%	0.3 to 1.5%	11903456
11895772	Fisherbrand Elite Stand, Accommodates 6 Pipettes					



# Fisherbrand® SureOne Pipet Tips

Fisherbrand™ SureOne™ pipet tips are a comprehensive line of universal fit pipet tips offered in many packaging configurations. Tips are available in volume range from 5 µl to 10 ml.

- Sterile product certified free of RNase, DNase, Pyrogen, Bioburden, PCR Inhibitors, Endotoxin and DNA
- Genomic tips and environmentally-friendly reload packaging also available
- Non sterile product certified free of RNase, DNase, DNA and Pyrogen
- Optimal fit with the Fisherbrand pipettes

## SureOne Filter Tips

Fisher Scientific Cat. No.	Tip Volume	Packaging	Color	Elite Pipette Fits Guide	
11917734 11997714	0.1 - 5 µl	Filtered Sterile Filter Sterile Reload	Clear	11815762	 Universal Fits, Micropoint Nano tip, graduated at 2.5 µl. Tip length 1.22 in, 3.10 cm
11903466 11907724	0.1 - 10 µl	Filtered Sterile Filter Sterile Reload	Clear	11815762 11825762 11835762 11815772 11855772	 Universal Fits, Micropoint tip, graduated at 2 µl. Tip length 1.23 in, 3.13 cm
11913466 11917724	0.1 - 10 µl	Filtered Sterile Filter Sterile Reload	Clear	11815762 11825762 11835762 11815772 11855772	 Universal Fits, Micropoint tip graduated at 2, 5, 10, 20 µl, Extended length. Tip length 1.80 in, 4.57 cm
11977714 11927724	0.1 - 20 µl	Filtered Sterile *Filter Sterile Reload	Clear	11815762 11825762 11835762 11815772 11855772	 Universal Fits, Micropoint tip graduated at 2, 5, 10, 20 µl, Extended length. Tip length 1.80 in, 4.57 cm
11943466 11937724	2 - 20 µl	Filtered Sterile **Filter Sterile Reload	Clear	11845762	 Universal Fits, Beveled tip. Tip length 2.01 in, 5.11 cm
11953466 11947724	10 - 100 µl	Filtered Sterile Filter Sterile Reload	Clear	11855762 11865762 11835772 11875772 11825772 11865772	 Universal Fits, Beveled tip. Tip length 2.01 in, 5.11 cm
11963466 11957724	20 - 200 µl	Filtered Sterile Filter Sterile Reload	Clear	11855762 11865762 11875762 11835772 11875772 11825772 11865772	 Universal Fits, Micropoint tip graduated at 10, 50, 100, 200 µl. Tip length 2.30 in, 5.84 cm
11987714 11967724	30 - 300 µl	Filtered Sterile Filter Sterile Reload	Clear	11855762 11865762 11875762 11835772 11875772 11825772 11865772 11845772 11885772	 Universal Fits, Micropoint tip graduated at 10, 50, 100, 200 µl. Tip length 2.30 in, 5.84 cm
11973466 11977724	100 - 1000 µl	Filtered Sterile Filter Sterile Reload	Clear	11885762	 Universal Fits, Micropoint tip graduated at 100, 200, 500, 1000 µl. Tip length 3.50 in, 8.89 cm
11927754	1 - 5 ml	Filtered Sterile, Bulk	Clear	NA	 Fits Gilson Pipetman and Gilson Style pipettes. Tip length 4.60 on, 11.70 cm

\*Use empty rack 10 - 20 µl: 02-707-445 \*\*Use empty rack 20 - 300 µl: 02-707-444

## SureOne Standard Tips

Fisher Scientific Cat. No.	Tip Volume	Packaging	Color	Elite Pipette Fits Guide	
11987724 11997724 11907734	0.1 - 5 µl	Bulk Racked Racked Sterile	Clear	11815762 11825762 11835762 11815772 11855772	 Universal Fits, Micropoint Nano tip, graduated at 2.5 µl. Tip length 1.22 in, 3.10 cm
11933416 11953416 10053014 11973416	0.1 - 10 µl	Bulk Racked Racked Sterile Reload	Clear	11815762 11825762 11835762 11815772 11855772	 Universal Fits, Micropoint tip, graduated at 2 µl. Tip length 1.23 in, 3.13 cm
11983416 11588402 10527014 11967714	0.1 - 10 µl	Bulk Racked Racked Sterile Reload	Clear	11815762 11825762 11835762 11815772 11855772	 Universal Fits, Micropoint tip graduated at 2, 5, 10, 20 µl, Extended length. Tip length 1.80 in, 4.57 cm
11933426 10678325 11963426 11578412 10124314 11983426 11903436 11913436	1 - 200 µl	Bulk Racked Racked Sterile Reload Bulk Racked Racked Sterile Reload	Clear Clear Clear Clear Yellow Yellow Yellow Yellow	11845762 11855762 11865762 11875762 11835772 11825772 11875772 11865772	 Universal Fits, Beveled tip. Tip length 2.01 in, 5.11 cm
11538422 11933436 11943436 11953436 11963436 11973436 11983436 11993436	1 - 200 µl	Bulk Racked Racked Sterile Reload Bulk Racked Racked Sterile Reload	Clear Clear Clear Clear Yellow Yellow Yellow Yellow	11845762 11855762 11865762 11875762 11835772 11825772 11875772 11865772	 Universal Fits, Thin Wall Micropoint tip. Tip length 1.97 in, 5.00 cm
11943446 11953446 11963446 12922521 11903446 11913446 11923446 11933446	1 - 200 µl	Bulk Racked Racked Sterile Reload Bulk Racked Racked Sterile Reload	Clear Clear Clear Clear Yellow Yellow Yellow Yellow	11845762 11855762 11865762 11875762 11835772 11825772 11875772 11865772	 Universal Fits, Beveled tip grated at 10, 50, 100 µl. Tip length 1.98 in, 5.04 cm
11993446 11903456 10003414 11538442	5 - 300 µl	Bulk Racked Racked Sterile Reload	Clear	11845762 11855762 11865762 11875762 11835772 11885772	 Universal Fits, Micropoint tip graduated at 10, 50, 100, 200 µl. Tip length 2.30 in, 5.84 cm
11548442 11568442 10164694 11588442 10778535 10537014 10492725 11963456	100 - 1250 µl	Bulk Racked Racked Sterile Reload Bulk Racked Racked Sterile Reload	Clear Clear Clear Clear Blue Blue Blue Blue	11885762	 Universal Fits, Micropoint tip graduated at 100, 200, 500, 1000 µl. Tip length 3.50 in, 8.89 cm
11987744	250 µl - 2.5 ml	Bulk	Clear	NA	 Fits with Rainin EDP2. Tip length 4.14 in, 10.53 cm
11997744	500 µl - 2.5 ml	Bulk	Clear	NA	 Fits Eppendorf and Biohit Style pipettes. Tip length 4.52 in, 11.48 cm
11648138	1 - 5 ml	Bulk	Clear	NA	 Fits Eppendorf and Biohit. Tip length 5.92 in, 15.04 cm
11937754	1 - 5 ml	Bulk	Clear	11895762	 Fits Fisherbrand Elite, Fisherbrand II and Thermo Scientific Finnpiquette Style. Tip length 5.80 in, 14.74 cm
11917754	1 - 5 ml	Bulk	Clear	NA	 Fits Gilson Pipetman and Gilson Style pipettes. Tip length 4.60 in, 11.70 cm
11947754	1-10 ml	Bulk	Clear	NA	 Fits Gilson Pipetman and Gilson Style pipettes. Tip length 6.32 in, 16.07 cm
11957754	1-10 ml	Bulk	Clear	11805772	 Fits Fisherbrand Elite, Fisherbrand II and Thermo Scientific Finnpiquette Style. Tip length 6.00 in, 15.24 cm

## SureOne Specialty Tips

Fisher Scientific Cat. No.	Tip Volume	Packaging	Color	Elite Pipette Fits Guide	
11927734 11937734 11367801	1 - 200 µl	Bulk Racked Racked Sterile	Clear	11845762 11855762 11865762 11875762	 Universal Fits, Gel Loading Tip, Outer Diameter 0.58 mm. Tip length 2.70 in, 6.86 cm
12992571	2 - 20 µl	Filtered Sterile	Clear	11845762	 Universal Fits, Gel Loading Tip, Outer Diameter 0.58 mm. Tip length 2.70 in, 6.86 cm
11967734	10 - 100 µl	Filtered Sterile	Clear	11855762 11865762	 Universal Fits, Gel Loading Tip, Outer Diameter 0.58 mm. Tip length 2.70 in, 6.86 cm
11977734 11987734 11997734	1 - 200 µl	Bulk Racked Racked Sterile	Clear	11845762 11855762 11865762 11875762	 Universal Fits, Extended Length Tip, 9 cm long. Tip length 3.57 in, 9.09 cm
12942581	20 - 200 µl	Filtered Sterile	Clear	11845762 11855762 11865762 11875762	 Universal Fits, Extended Length Tip, 9 cm long. Tip length 3.57 in, 9.09 cm
11917744 11927744 11937744	1 - 200 µl	Bulk Racked Racked Sterile	Clear	11845762 11855762 11865762 11875762 11835772 11825772 11875772 11865772	 Universal Fits, Low Retention, Large Orifice Tip, Inner Diameter 1.50 mm. Tip length 1.96 in, 4.98 cm
11947744	20 - 200 µl	Filtered Sterile	Clear	11845762 11855762 11865762 11875762 11835772 11825772 11875772 11865772	 Universal Fits, Low Retention, Large Orifice Tip, Inner Diameter 1.50 mm. Tip length 1.96 in, 4.98 cm
12992581 11957744 11967744	100 - 1000 µl	Bulk Racked Racked Sterile	Clear	11885762	 Universal Fits, Low Retention, Large Orifice Tip Graduated 100, 200, 500, 1000 µl, Inner Diameter 1.50mm. Tip length 3.50 in, 8.89 cm
11977744	100 - 1000 µl	Filtered Sterile	Clear	11885762	 Universal Fits, Low Retention, Large Orifice Tip, Graduated 100, 200, 500, 1000 µl, Inner Diameter 1.50 mm. Tip length 3.50 in, 8.89 cm

## SureOne Reload Empty Racks

Fisher Scientific Cat. No.	Tip Volume	Packaging	Color	Elite Pipettes Fits Guide	
11973456 11983456 11993456	10 - 20 µl 20 - 300 µl 1250 µl	Empty Rack Empty Rack Empty Rack		NA NA NA	Inserts not included



# Fisherbrand® SureOne Pipet Tips







## SureOne Filter Tip Reloading Steps



**Step 1** While grasping below the lowest insert, lift sleeve up off rack and transfer entire sleeve to an empty rack. Make sure to align sleeve knockouts with latch and rear lid posts.

**Step 2** Push down with even pressure on the middle of the top insert, which will slide bottom insert into empty tray. Plastic sleeve is not intended to rest on any counter top without a rack underneath it.

## SureOne Pipet Tips Packaging Options

						
	<b>Bulk</b>	<b>Racked</b>	<b>Racked Sterile</b>	<b>Reload</b>	<b>Filter Tip Reload</b>	<b>Empty Rack</b>
Pack Format	1000 tip x bag	96 tips x tray 10 trays x pk	96 tips x tray 10 trays x pk	96 tips x tray 10 trays x pk	*2 Stacks of 480 tips 96 tips x tray 10 trays x pk	Inserts not included 10 empty trays x pk
Case Format	10 bags x cs	5 pks x cs	5 pks x cs	5 pks x cs	5 pks x cs	5 pks x cs

\*1 ml Filter Tip Reload: 56 tips x tray, 10 trays x pk, 5 pks x cs

© 2016 Thermo Fisher Scientific Inc. All rights reserved.  
Trademarks used are owned as indicated at [www.fishersci.com/trademarks](http://www.fishersci.com/trademarks).

**Austria:** +43(0)800-20 88 40 **Belgium:** +32 (0)56 260 260 **Denmark:** +45 70 27 99 20  
**Germany:** +49 (0)2304 9325 **Ireland:** +353 (0)1 885 5854 **Italy:** +39 02 950 59 478  
**Finland:** +358 (0)9 8027 6280 **France:** +33 (0)3 88 67 14 14 **Netherlands:** +31 (0)20 487 70 00  
**Norway:** +47 22 95 59 59 **Portugal:** +351 21 425 33 50 **Spain:** +34 902 239 303  
**Sweden:** +46 31 352 32 00 **Switzerland:** +41 (0)56 618 41 11 **UK:** +44 (0)1509 555 500

 **Fisher Scientific**  
A Thermo Fisher Scientific Brand