


## Pipettes / Pipette fillers

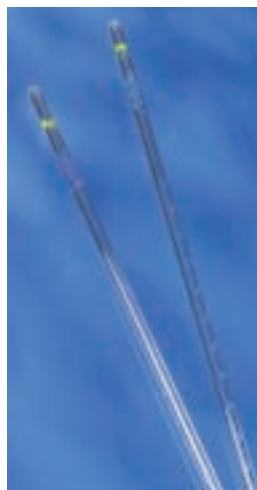

**Pipettes,**  
soda lime glass,  
graduated, BS700, ISO835,  
Class A

**Fisherbrand®**

- Blue ceramic graduations
- Supplied in packs of 2

Catalogue No	Cap., mL	Grads., mL	Tol., mL
Type 1 Calibrated to deliver from zero to any graduation line			
FB50231	1	0.01	±0.01
FB50233	2	0.02	±0.02
FB50235	5	0.05	±0.05
FB50237	10	0.01	±0.10
FB50239	25	0.02	±0.20
Type 2 Calibrated to deliver from any graduation line to top jet			
FB50241	1	0.01	±0.01
FB50243	2	0.02	±0.02
FB50245	5	0.05	±0.05
FB50247	10	0.01	±0.10
FB50249	25	0.02	±0.20



**Pipettes, glass, Pasteur**

**Fisherbrand®**

- Supplied in packs of 1,000 (4 x 250) without teats

Catalogue No	Nominal length, mm
FB50251	150
FB50253	230
Catalogue No	Description
PMK-302-500E	Pipette teats, PVC, pack of 100

**Pipette filler, Bibbette**

- PTFE and silicone rubber valve assembly controls the vacuum created by simply squeezing the unique squared rubber bulb
- Pipettes from 0.1 to 100mL can be easily and accurately filled and emptied
- Blow out pipettes are catered for by the depression of a small thumb operated squeeze ball
- Replaceable PTFE membrane filter ensures liquid is not accidentally drawn into the unit
- Filter and holding cone are autoclavable
- Supplied with two membrane filters, full instructions and six month warranty



Catalogue No  
PMR-600-010Y

Spare Part	Description
Catalogue No	
PMR-600-030C	Spare membrane filters, pack of 10



**Pipette, bulbs,**  
rubber teats

- Supplied in packs of 20

Catalogue No	Neck bore, mm	Nominal cap., mL
PMR-210-011Y	5	1
PMR-210-031C	5	1.5
PMR-210-051T	6	2.5
PMR-210-071N	6	3.75
PMR-220-G	Mixed pack of 5 of each of the above	

**Fisherbrand®**

**Pipette fillers**



- Moulded rubber bulb, approximately 60mL capacity, with three glass ball valves which can be controlled by finger pressure
- Can be used for all liquids except those which attack rubber
- Can be used on bulb or cylindrical pipettes from 2 to 50mL
- Do not allow liquids to come into contact with the valves or inside of the bulbs, this may impair performance

Catalogue No	Colour
FB51420	Black
FB51421	Blue
FB51422	Red
FB51423	Yellow

**Pipette filler, diji™**

**NEW**

- Designed for handling laboratory liquids
- Accuracy enhanced due to smooth handling of head
- Wrist and finger tension eliminated
- More secure - better grip than standard pipette fillers
- Easy to attach, clean and change pipettes
- Teflon® sheathing eliminates broken glass and stabbing injuries commonly caused by breakage of pipettes
- Fluid cannot be accidentally sucked into the mechanism



Catalogue No	Volume, mL
PMR-300-250M	10
PMR-300-100Y	25